

# ERBE



Manufacturer:

**ERBE Elektromedizin GmbH**  
Waldhörnlestr. 17  
72072 Tübingen, Germany  
P.O. Box 14 20  
72004 Tübingen, Germany  
Tel.: (+49) (70 71) 755-0  
Fax: (+49) (70 71) 755-179  
Email: support@erbe-med.de

Distributed by (USA Territory):

**ERBE USA, Inc.**  
2225 Northwest Parkway  
Marietta, GA 30067  
Customer Support 800-778-ERBE (3723)  
Tel.: (770) 955-4400  
Email: sales@erbe-usa.com

07.10




Printed in Germany  
© ERBE Elektromedizin GmbH

20139-174	20150-220	20150-238
20139-178	20150-225	20150-239
20139-181	20150-226	20150-241
20139-185	20150-230	
20139-186	20150-231	

30150-001

**Single use only!** Sterility not guaranteed if package has been opened or damaged. Discard if package is damaged.

**STERILE EO**

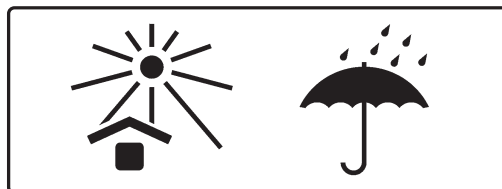

**USA**

## Applicators and Probe

### IMPORTANT!

Please read all information carefully.

These notes on use do not replace the user manual of the Water Jet System used! Read the user manual of the Water Jet System and ask ERBE or your distributor in case of doubt!

**Caution!** Federal law (USA) restricts this device to sale by or on the order of a physician!

### Purpose

The applicators and probe are used with designated Water Jet Units for the cutting and dissection of soft tissue in clinical applications as outline in the intended use of the Water jet Units.

### Intended Use

The applicators and probe are intended for use with the following designated ERBE Equipment. For further information see User Manual of the specific ERBE equipment.

- **Applicators for ERBE Helix Hydro-Jet:** Applicators with Art. No. 20139-xxx. (**Note:** Only ERBE Fluid cartridge Art. No. 20139-083 is to be used in conjunction with these applicators).
- **Applicators and Probe for ERBEJET 2:** Accessories with Art. No. 20150-xxx. (**Note:** Only ERBE Pump cartridge Art.

No. 20150-300 is to be used in conjunction with these accessories).

Applicators for the ERBE Helix Hydro-Jet that have a suction hose (with suction) can be connected to the suction of the Helix Hydro-Jet.

Applicators for the ERBEJET 2 that have a suction hose (with suction) can be connected to the suction module, Model ESM 2 (**Note:** Alternate suction devices can be used but the suction (vacuum) pressure should not exceed -11.6 psi (-800 mbar)).

## Safety Instructions

### WARNING!

This product may only be used by medical personnel who have been trained per these instructions.

Before use, check that the product is intact. Verify that there are no signs of physical damage.

If damaged, do not use this product!

When using and handling the product, the involved medical staff should wear personal protective equipment (e.g. gown, faceshield or eye protection, mask, gloves, etc.) to minimize the risk of exposure to infectious agents due to undesirable/unexpected spraying or splashing, mishandling of the contaminated accessory, etc.

When using an applicator with suction, maintain adequate and constant suction as necessary as well as monitor the suction tubing to minimize the possibility of blockage.

Never lay this product on the patient or in his/her direct vicinity.

Protect this product from any form of mechanical damage! Do not throw! Do not use force!

ERBE recommends that the product not be reprocessed after use. High temperatures in particular (e.g. in the autoclave) can affect material properties and product function making its Intended Use no longer possible.

ERBE Elektromedizin expressly warns against modifying the product. Any modification exempts ERBE Elektromedizin from any and all liability.

## How to Use

### Check the expiration date

An expiration date is printed on the packaging of the product. Do not use any product beyond its expiration date!

### Control and open packaging

Use this product only if the packaging has not been opened or damaged. If the packaging is open or damaged, the product is no longer sterile and must not be used. In this case, throw the product away!

Use aseptic techniques when removing the product from the packaging!

### Connecting the product

To connect the applicators or probe, please read the relevant section in the following designated ERBE equipment User Manual.

- For applicators with Art. No. 20139-xxx: ERBE Helix Hydro-Jet User Manual
- For applicators and probe with Art. No. 20150-xxx: ERBEJET 2 User Manual (Note: The probe connects and is used like an applicator.)